SIMPLETaste

BATHROOM SCALE

1byone Products Inc.
1230 E Belmont Street, Ontario, CA 91761
Customer Service: +1 909-391-3888
www.1byone.com


Instruction Manual/Bedienungsanleitung/Manuel d'instructions
Manual de instrucciones/Manuale di istruzioni

Introduction

Thank you for purchasing SimpleTaste Bathroom Scale. This instruction manual contains important information about safety, usage, and disposal. Use the product as described and keep this manual for future reference. If you sell this scale or pass it on, also give this manual to the new owner.

Appliance Description

1. LCD display
2. Battery compartment
3. Unit button

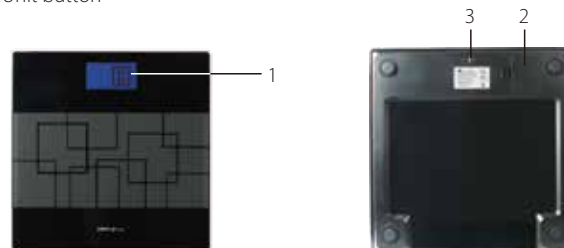Safety Instructions

Intended Use

This product is only authorized for domestic use and not for commercial or industrial purposes. It is only intended for applications described in these instructions.

Safety Information

- Do not immerse the scale in water or other liquids.
- Do not place the scale on carpets or soft surfaces. Only place on solid and flat surfaces.
- Do not step on the scale with wet feet and do not use the scale if the platform is moist. Moisture will cause the scale's surface to become

slippery.

- When using the scale for the first time, or when changing its location, it should be turned on, off, and on again. It's surface should also be free of any objects. Failure to do this may result in inaccurate measurements.
- Always keep batteries out of the reach of infants and young children to prevent them from being swallowed. If swallowed, take immediate action and contact medical authorities for help.

Cleaning

Clean the scale with a slightly damp cloth. Do not use any chemical/abrasive cleaning agents.

Storage

If you do not intend to use the scale for a long period of time, remove the batteries and store the scale in a clean, dry place away from direct sunlight.

Operation

Included Items

- 1 x Bathroom scale
- 3 x AAA batteries


Before using the scale for the first time, check the package contents to ensure all parts are free of visible damage.

Inserting batteries


To insert the batteries, open the battery compartment on the rear of the scale and remove the compartment cover. Insert 3 AAA batteries, paying attention to the correct polarity.

Select unit before use

1. Wait until "0.0lb" appears.

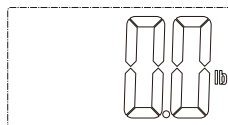
2. Press the UNIT button on the rear of the scale to select the desired unit (kg, lb, or st).

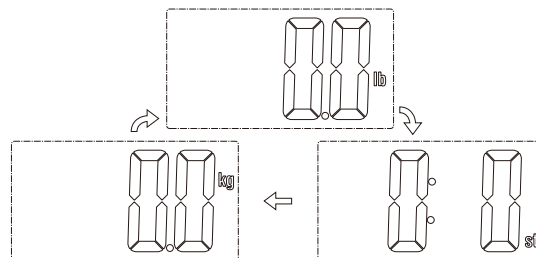
Note: After selecting the unit, place the scale on a solid and flat surface and start the measurement as described below. The scale will shut off automatically after 10 seconds of inactivity.

Measuring (with lb unit selected)

1. Place the scale on a solid and flat surface.
2. Stand on the scale to turn it on.

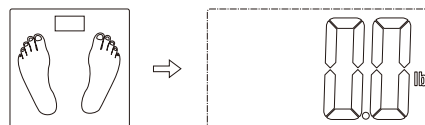
3. Stand still and the display will blink three times presenting the confirmed weight measurement. The scale will then turn off after 10 seconds.

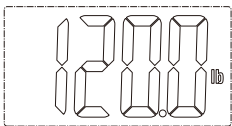
Note: If you want to measure again, please step off completely and then step back on the scale.

Indication messages

- "Lo" is shown when battery power is low. The scale will then turn off after a few seconds. Please replace with new batteries.

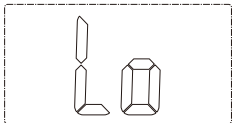- "Err" is shown when the scale is overloaded. Please remove the weight on the scale immediately to protect the sensors.

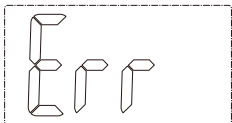Specifications

Capacity: maximum 180kg/400lb

Graduation: 0.05kg/0.1lb

Weight unit: kg/lb/st:lb

Minimal weight: 5kg

Batteries: 3 x AAA

Disposal

Disposal of the appliance

Under no circumstances should you dispose of the appliance in normal domestic waste. This product is subject to the provisions of European Directive 2012/19/EU.

Dispose of the appliance via an approved disposal company or your municipal waste facility. Please observe currently applicable regulations. Please contact your waste disposal center if you need any further information.

Disposal of batteries

Used batteries may not be disposed of in household waste.

All consumers are statutorily required to dispose of batteries at a collection point provided by their local municipality or retail store.

The purpose of this obligation is to ensure batteries are disposed of in a non-polluting manner. Only dispose of batteries when they are fully discharged.

The appliance's packaging is made from environmentally friendly material and can be disposed of at your local recycling plant.

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